

User Guide

Doherty Paper Roll Holder - CHA02

DOHERTY

Thank you for purchasing a paper roll holder for your Doherty treatment chair. The paper roll holder has been designed specifically for:

CHE01 - Doherty Phlebotomy chair

CHE02 – Doherty Variable Height Treatment Chair


Kit of parts:

- 1x Paper roll bracket
- 1x Paper roll Bar
- 2x Ball end caps
- 2x M6 Button head bolts
- 1x Allen key

1. Offer the paper roll holder bracket up to the rivnut holes on the backrest.
2. Once lined up use the M6 bolts to hold it in place, using the Allen key to tighten up and secure.
3. To attach a paper roll, remove one of the ball end caps and slide the bar out from the bracket. Locate the roll and then re attach the ball end cap.



Before use always ensure both the bolts and ball end caps are fully located.

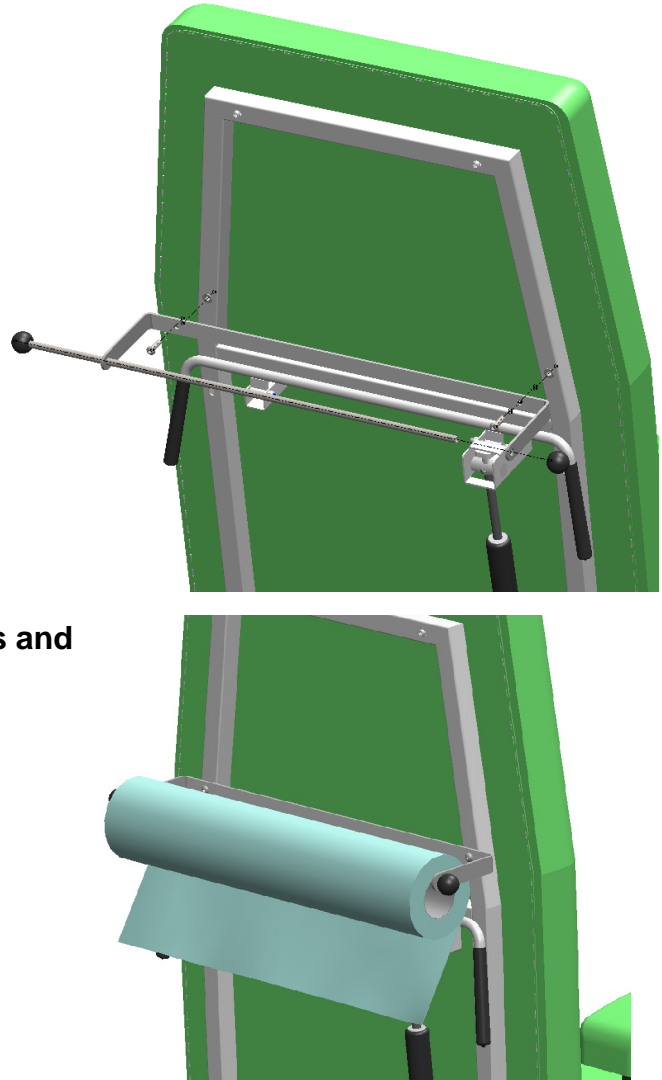

This user guide should be read in conjunction with the Doherty Treatment Chair user guide. Please refer to the Chair user guide for cleaning procedures, maintenance and warranty information.

A copy of the Doherty Treatment Chair user guide can be obtained from a member of the Sidhil customer services team on the telephone number shown below.

Sidhil Limited Sidhil Business Park Holmfield Halifax West Yorkshire HX2 9TN

Tel: +44 (0)1422 233030 Fax: +44 (0)1422 233040 Email: info@sidhil.com
Sales Tel: +44 (0)1422 233030 Sales Fax: +44 (0)1422 233040

A Member of the Siddall & Hilton Ltd Group of Companies

SIDHIL

DOHERTY

This document was created with Win2PDF available at <http://www.daneprairie.com>.
The unregistered version of Win2PDF is for evaluation or non-commercial use only.